




## 维生素E异构体分析

AN0113

色谱柱: SenSil™ C30, 5 $\mu$ m, 4.6 x 250mm  
规格: 250 x 4.6 mm ID  
货号: 2061055.s2546  
洗脱方式: Isocratic  
流动相: Methanol/Water (99:1)  
流速: 1.0mL/min  
柱温: 15°C  
进样体积: 1 $\mu$ L  
检测器: UV@295nm

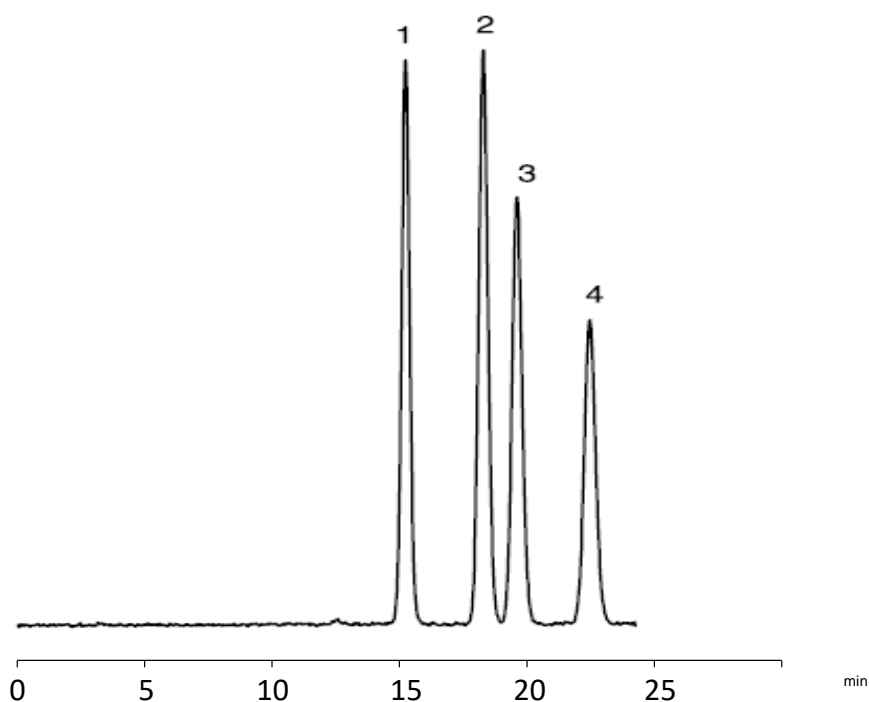

### 化合物名称:

1.  $\delta$ -生育酚
2.  $\gamma$ -生育酚
3.  $\beta$ -生育酚
4.  $\alpha$ -生育酚